



• **GenoXtract[®] fleXT**

Extraction and PCR set-up for FluoroType[®] assays* directly from original patient specimens

- Walk-away extraction automation: Minimal and easy workflow set-up to maximize your hands-off time
- Accurate and efficient: Parallel and highly reliable processing of all samples under continuous deviation monitoring control
- Safe and clean: Loading outside the pipetting area to avoid contamination and allow ergonomic interaction

GenoXtract[®] fleXT is an automated extraction device, with a flexible option to set up diagnostic PCR. Up to 96 samples can be directly extracted from patient specimens for IVD purposes.

*Please contact your local representative for the current list of validated assays

Minimum hands-on, maximum walk-away time


GenoXtract® fleXT is your reliable, automated extraction system directly from the original patient sample with no need for extensive preparation prior to extraction.

Barcode-supported robotic sorting significantly reduces manual intervention. Optionally combine extraction with PCR set-up, for a convenient walk-away solution tailored to the dynamic laboratory environment.

Active liquid sensing for continuous volume monitoring with external loading and simplified data transfer to FluoroCycler® XT ensures sample extraction is more reliable and accurate than ever.

- ▶ Seamless workflow in combination with FluoroCycler® XT
- ▶ Minimal manual pipetting before and after extraction, frees up time for other tasks
- ▶ Peace of mind: accurate results supported by integrated deviation control

From patient sample to data analysis

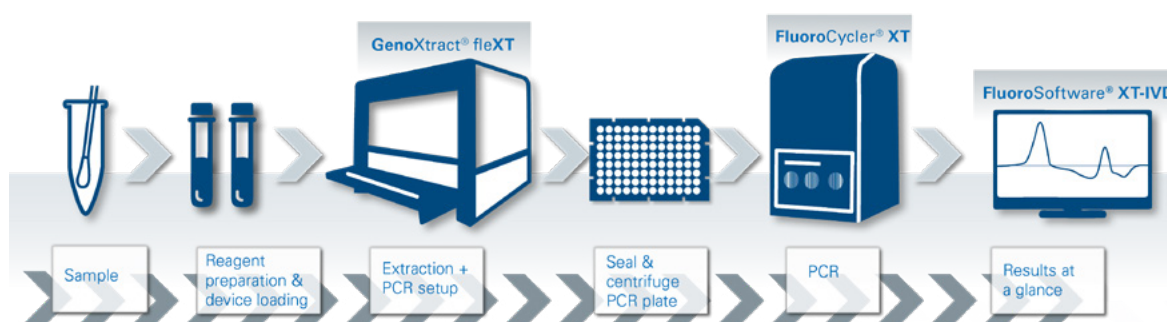

Specifications

Number of samples:	up to 96 with 8 pipetting channels
Pipetting method:	Micro Piston Pump (MPP) with Liquid Level Detection
Pipetting volume:	10–1000 µl
Contamination control:	integrated UV lamp
Dimensions:	112.4 × 90.3 × 79.5 cm (W × H × D)
Weight:	approx. 140 kg

Please contact your local representative for availability in your country. Not for sale in the USA.
FluoroType®, FluoroCycler® and GenoXtract® are registered trademarks of Bruker corporate group
GenoXtract® fleXT Hamilton Bonaduz AG, Via Crusch 8, CH-7402 Bonaduz, Switzerland
FluoroCycler® XT Hain Lifescience, Hardwiesenstrasse 1, 72147 Nehren, Germany

Issued by

● **Hain Lifescience GmbH**

Hardwiesenstrasse 1
72147 Nehren - Germany
Phone +49 (0) 7473-9451-0

info.mdx.de@bruker.com - www.hain-lifescience.de