



## • **GenoXtract<sup>®</sup> fleXT**

### **Extraction and PCR set-up for FluoroType<sup>®</sup> SARS-CoV-2 assays\* directly from original patient specimens**

- Walk-away extraction automation: Minimal and easy workflow set-up to maximize your hands-off time
- Accurate and efficient: Parallel and highly reliable processing of all samples under continuous deviation monitoring control
- Safe and clean: Loading outside the pipetting area to avoid contamination and allow ergonomic interaction

GenoXtract<sup>®</sup> fleXT is an automated extraction device, with a flexible option to set up diagnostic PCR. Up to 96 samples can be directly extracted from patient specimens for IVD purposes.

\*Please contact your local representative for the current list of validated assays


# More samples, more power against COVID-19

GenoXtract® fleXT is your reliable, automated extraction system for SARS-CoV-2 diagnostics\* directly from the original patient sample with no need for extensive preparation prior to extraction.

Barcode-supported robotic sorting significantly reduces manual intervention. Optionally combine extraction with PCR set-up, for a convenient walk-away solution tailored to the dynamic laboratory needs during and beyond the pandemic.

Active liquid sensing for continuous volume monitoring with external loading and simplified data transfer to FluoroCycler® XT ensures sample extraction is more reliable and accurate than ever.

- ▶ Seamless workflow across SARS-CoV-2 portfolio in combination with FluoroCycler® XT
- ▶ Freeing up time for other tasks by minimal manual pipetting before and after extraction
- ▶ Peace of mind: accurate results supported by integrated deviation control

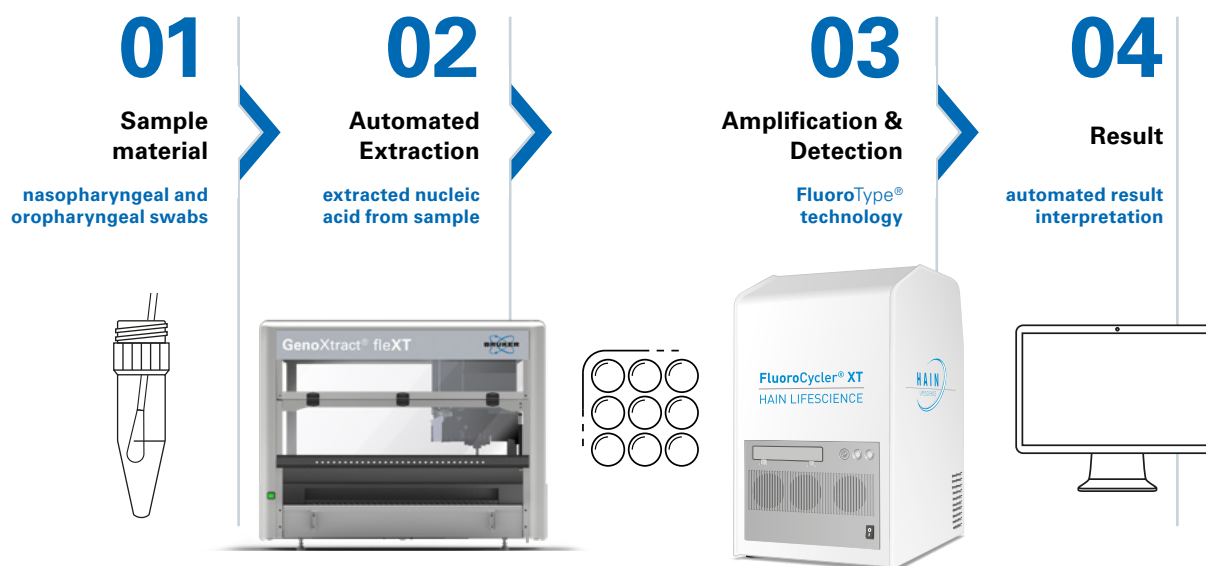

## Specifications

<b>Order number:</b>	1883774
<b>Number of samples:</b>	up to 96 with 8 pipetting channels
<b>Pipetting method:</b>	Micro Piston Pump (MPP) with Liquid Level Detection
<b>Pipetting volume:</b>	10–1000 µl
<b>Contamination control:</b>	integrated UV lamp
<b>Dimensions:</b>	112.4 × 90.3 × 79.5 cm (W × H × D)
<b>Weight approx:</b>	140 kg

Please contact your local representative for availability in your country. Not for sale in the USA. FluoroType®, FluoroCycler® and GenoXtract® are registered trademarks of Bruker corporate group.



Hardwiesenstrasse 1  
72147 Nehren - Germany  
Phone +49 (0) 7473-9451-0



[info.mdx.de@bruker.com](mailto:info.mdx.de@bruker.com) - [www.hain-lifescience.de](http://www.hain-lifescience.de)