

IVD



## FluoroType<sup>®</sup> SARS-CoV-2/Flu/RSV

- Real-Time PCR for differentiated detection of SARS-CoV-2, Influenza A virus, Influenza B virus and RSV. All in one test.

### Diagnostics for the Winter Needs

The diagnostic needs for the regular flu season during times of COVID-19 pose new challenges to diagnostic laboratories and doctors. With very similar symptoms of viral respiratory infection, a reliable test as FluoroType<sup>®</sup>SARS-CoV-2/Flu/RSV allows rapid differentiation.

## Easy detection of viruses of respiratory indications

- Multiplex Real-Time PCR
- 5 analytes simultaneously
- Differentiates between Flu A and Flu B
- Directly from respiratory samples (nasopharyngeal/ oropharyngeal swabs) and saliva
- Results in less than two hours
- High test capacity

### Targets

#### Flu A

Subtypes H1N1 and H3N2

#### Flu B

Lineages Victoria and Yamagata


#### RSV

RSV-A and RSV-B

#### SARS-CoV-2

#### Internal Control

### Result Example

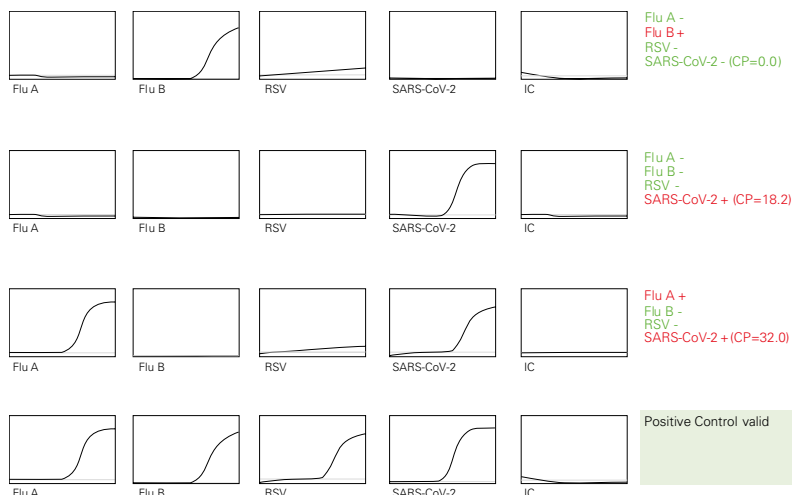

### Results at a glance

Automated RT-PCR on FluoroCycler® XT in combination with the fully automated test analysis gives the result at a glance, easy to use, easy to interpret.

The FluoroType® SARS-CoV-2/Flu/RSV kit is combined with an encapsulated RNA Internal Control (Universal Internal Control 2), which mimics lysis of viral particles and serves as an extraction, reverse transcription and amplification control.



### Detection of common winter viruses: Flu, RSV and SARS-CoV-2

- Cover your winter related respiratory viruses in only one test
- Improve your diagnostic time to result during the Flu season
- Distinguish your positive results at one glance
- Sample material: swab or saliva

### Order Information:

FluoroType® SARS-CoV-2/Flu/RSV  
**Part No. 1880015**

Universal Internal Control 2  
**Part No. 402960**

FluoroCycler® XT  
**Part No. 2017**

GenoXtract®  
**Part No. 8.31.01**

GXT NA Extraction Kit  
**Part No. 12.08.02**

Please contact your local representative for availability in your country. Not for sale in the USA.  
FluoroType®, FluoroCycler® and GenoXtract® are registered trademarks of Hain Lifescience GmbH – a Bruker company



● **Hain Lifescience GmbH**

Nehren · Germany  
Phone +49 (0) 7473-9451-0

[info.mdx.de@bruker.com](mailto:info.mdx.de@bruker.com)  
[www.hain-lifescience.de](http://www.hain-lifescience.de)